

## 18.0.Extra Oral

### 18.1. Multi Adjustable Facemasks



The adaptability of the Adjustable Protraction Face Mask allows for directional force application by adjusting the elastics engagement bar to deliver protractile force in the desired plane. Fully adjustable forehead and chin rests allow for adaptation to most patient's anatomical structure. The angle of traction is then determined by adjusting the elastics engagement bar to deliver horizontal, upward or downward forces as necessary. Face mask adjustments are made with an Allen wrench (included). Padded forehead and chin rests offer comfort for extended wear and may be cleaned. Traction force application recommends 1/2" elastics with 16 oz. pull ratings, but alternate elastics may be used if differing forces are desired. Sold individually;

replacement parts available separately. Include

**AE320-001** - Multi Adjustable Face Mask Dynamic Large

**AE320-002** - Multi Adjustable Face Mask Dynamic Medium

**AE320-003** - Multi Adjustable Face Mask Dynamic Small

### 18.2. Reverse Pull Face Mask


The Reverse Pull Face Mask is available in two sizes and offers comfortable anchorage for protraction. The stainless steel frame is softened by hypoallergenic chin and forehead pads to minimize irritation and maximize patient acceptance. Bilateral and longitudinal adjustability allow for patient customization. All adjustments are made with an Allen wrench (included with each mask). Each mask also includes: 1 forehead pad, 2 chin pads and replacement screws. Sold individually; replacement parts available separately.



**AE320-010** - Adjustable Reverse Pull Headgear Small

**AE320-011** - Adjustable Reverse Pull Face Mask Large

### 18.3. Standard Face bows

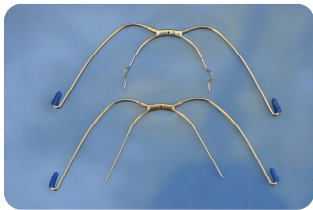

In bright stainless steel. Five sizes plus Universal - Available in short or medium outer bow sizes

#### Face bows

TYPE	SHORT	MEDIUM
<i>SS Face-bow Universal</i>	AE320-020	AE320-030
<i>SS Face-bow Sz 1</i>	AE320-021	AE320-031
<i>SS Face-bow Sz 2</i>	AE320-022	AE320-032
<i>SS Face-bow Sz 3</i>	AE320-023	AE320-033
<i>SS Face-bow Sz 4</i>	AE320-024	AE320-034
<i>SS Face-bow Sz 5</i>	AE320-025	AE320-035

## 18.4. Lip Bumpers



SIZE	CODE
<i>Lip Bumper Universal</i>	AE320-040
<i>Lip Bumper Sz.1</i>	AE320-041
<i>Lip Bumper Sz.2</i>	AE320-042
<i>Lip Bumper Sz.3</i>	AE320-043
<i>Lip Bumper Sz.4</i>	AE320-044
<i>Lip Bumper Sz.5</i>	AE320-045

Lip Bumpers are precision shaped to relieve lip and cheek pressures and allow teeth to unravel quickly. Bumpers feature colored lip pad and include small ligating hooks for securing in place. Sizes 1-5\* and Universal, .045" diameter wire - Two per pack.

## 18.5.High Pull Head caps



### Head caps

TYPE	SMALL	LARGE
<i>Black</i>	AE320-050	AE320-053
<i>Royal Blue</i>	AE320-051	AE320-054
<i>Red</i>	AE320-052	AE320-055

Two sizes of head caps, each with two placements loops, ensures accurate fit for patients. Pkg. 5.

## 18.6.Cervical Neck pads



TYPE	CODE
<i>Denim</i>	AE320-065
<i>Confetti</i>	AE320-066
<i>Pink</i>	AE320-067
<i>Royal Blu</i>	AE320-068

Two placement loops, for adjustability, ensure accuracy of fit. Manufactured with plush padding for comfort.

## 18.7. Safety Release Modules



Three force ranges to provide dependable retraction - Specially designed retention clip releases at 4-5 lbs. of outward force for patient safety. Safety Release Modules have been re-designed for greater strength and durability. **CODE** AE320-071