



23.2. Mini Mold

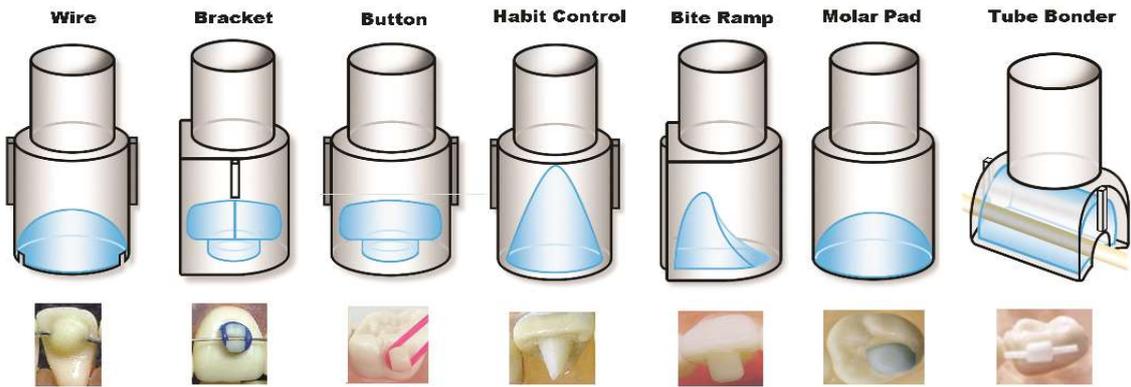

Light-cure your own resin attachments in seconds!

Make fixed retainers, tongue cribs, and bite-opening appliances without impressions. Mini Molds™ do the job quickly and easily with 12 different tips. Simply fill the tip, place on tooth and cure. In seconds, Mini Mold™ tips precisely cast a light-cured resin appliance. Our new Cuspid Bite Ramp is ideal for immediate bite opening when the overjet is in excess of 5mm which is a superior alternative when preventing any possible trauma to the incisors. These cured resin attachments are more comfortable on the tongue and very easy to equilibrate. Please visit our web site for Instruction

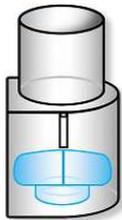


AE700-111



BRACKET BONDER

Forms a single wing bracket for simple correction. Splinting loose teeth due to traumatic injury.

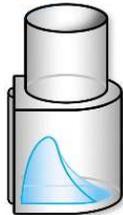


AE700-117 # 3 mm.
AE700-117L # 5 mm.



BITE RAMP

Forms a bite opening ramp on lingual of the upper incisors. Use to adjust the occlusion as needed.

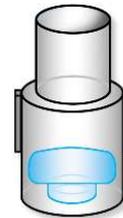

AE700-104



BUTTON BONDER

Forms a lingual button without impressions. Can be placed anywhere a button is needed.



AE700-116



HABIT CONTROL

Forms a small conical spike on the lingual of the incisors. Corrects thumb sucking and tongue thrusting with minimal chair time.




AE700-107



MOLAR PAD

Forms a precise spherical mound on the posterior cusps. Compliments the anterior bite-ramp for opening the bite from the posterior teeth.

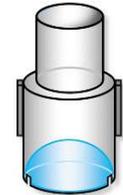

Small AE700-109
Large AE700-101



WIRE BONDER

Forms a Precise spherical resin mound over an arch wire. Ideal for making direct retainers in minutes. Excellent for rotational corrections



AE700-123



TUBE BONDER

Forms a tube for easy placement on erupting second molars. Use on the gingival of molar to hold a retainer clasp from slipping. Temporary tubes for relapse or simple cases

AE700-100 MINI MOLD STARTER KIT

Kit Contains:
Light-Cure Your Own Resin Attachments In Seconds! Make Invisible Aesthetic Attachments in seconds. Mini-Molds do the job quickly and easily. Just fill the mold (choose from 7 different treatment aids), place on tooth, and cure. In seconds, Mini-Mold tips precisely cast a light-cured resin appliance. Components Habit Control (3) Anterior Bite-Ramp 5mm (3) Anterior Bite-Ramp 3mm (3) Posterior Molar Pad (3) Wire Bonder small (3) Wire Bonder large (3) Button Bonder (3) Bracket Bonder (3) Tube Bonder (3) Bonder handle (2).

